
Neu ab:
New since:
Nouveau depuis:
Nuevo desde:

11.2004

sirona
The Dental Company

ORTHOPHOS Plus / Plus Ceph

Einbauanleitung für Ausbausatz Transversal-Schnitte
Installation Instructions for Transverse Section Attachment Kit
Instructions de montage pour jeu d'extension coupes transversales
Instrucciones de montaje para juego de ampl. para cortes transversales

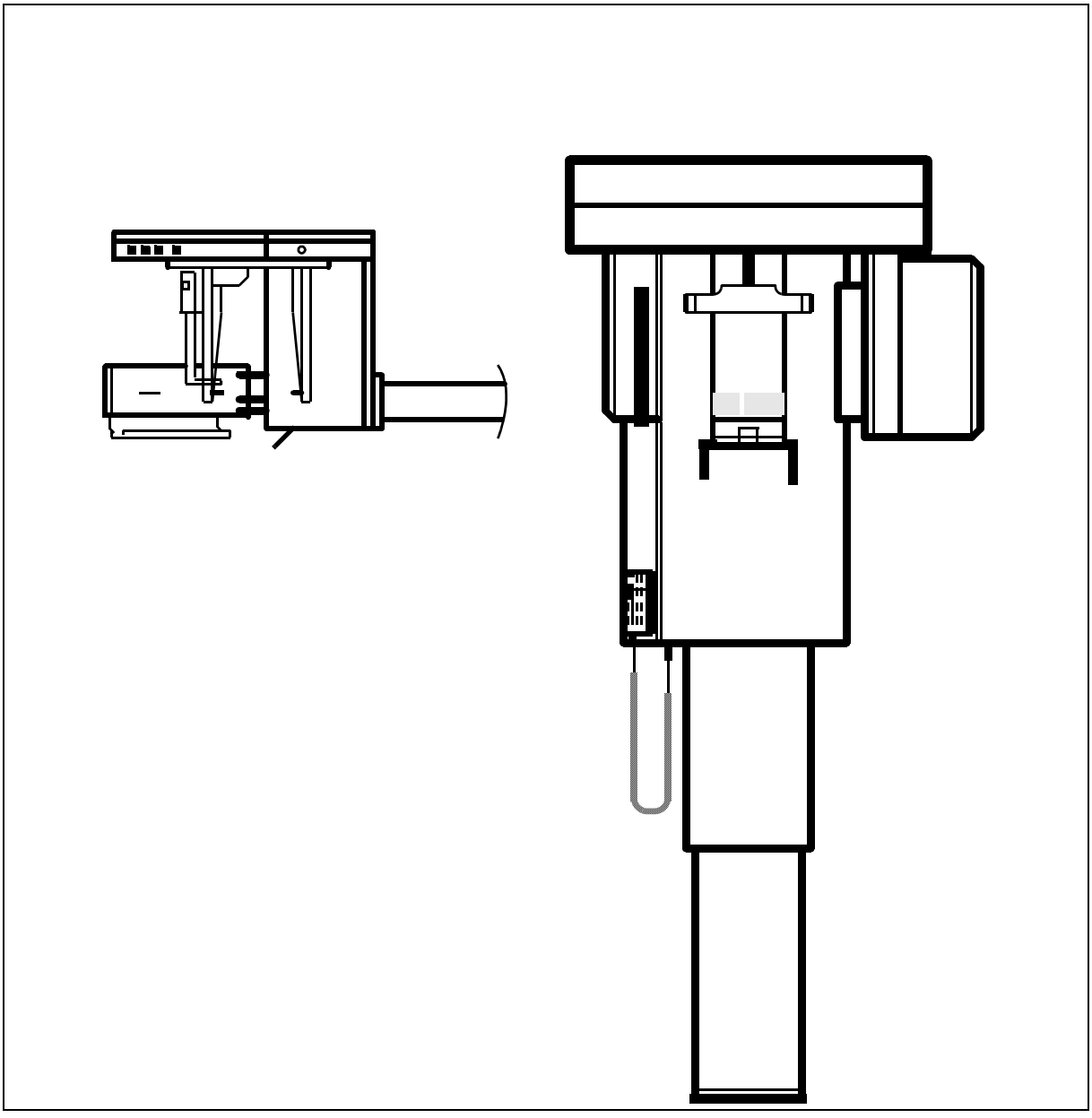ACHTUNG

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Sirona Dental Systems GmbH

Fabrikstraße 31
D-64625 Bensheim
Germany
www.sirona.com

in the USA:
Sirona Dental Systems LLC
4835 Sirona Drive, Suite 100
Charlotte, NC 28273
USA

in Canada:
Sirona Canada
3250 Ridgeway Drive - Unit 5
Mississauga, Ontario L5L 5Y6
Canada


Bestell-Nr. **54 71 607 D3297**

New since: 11.2004

sirona
The Dental Company

ORTHOPHOS Plus / Plus Ceph

Installation Instructions for Transverse Section Attachment Kit

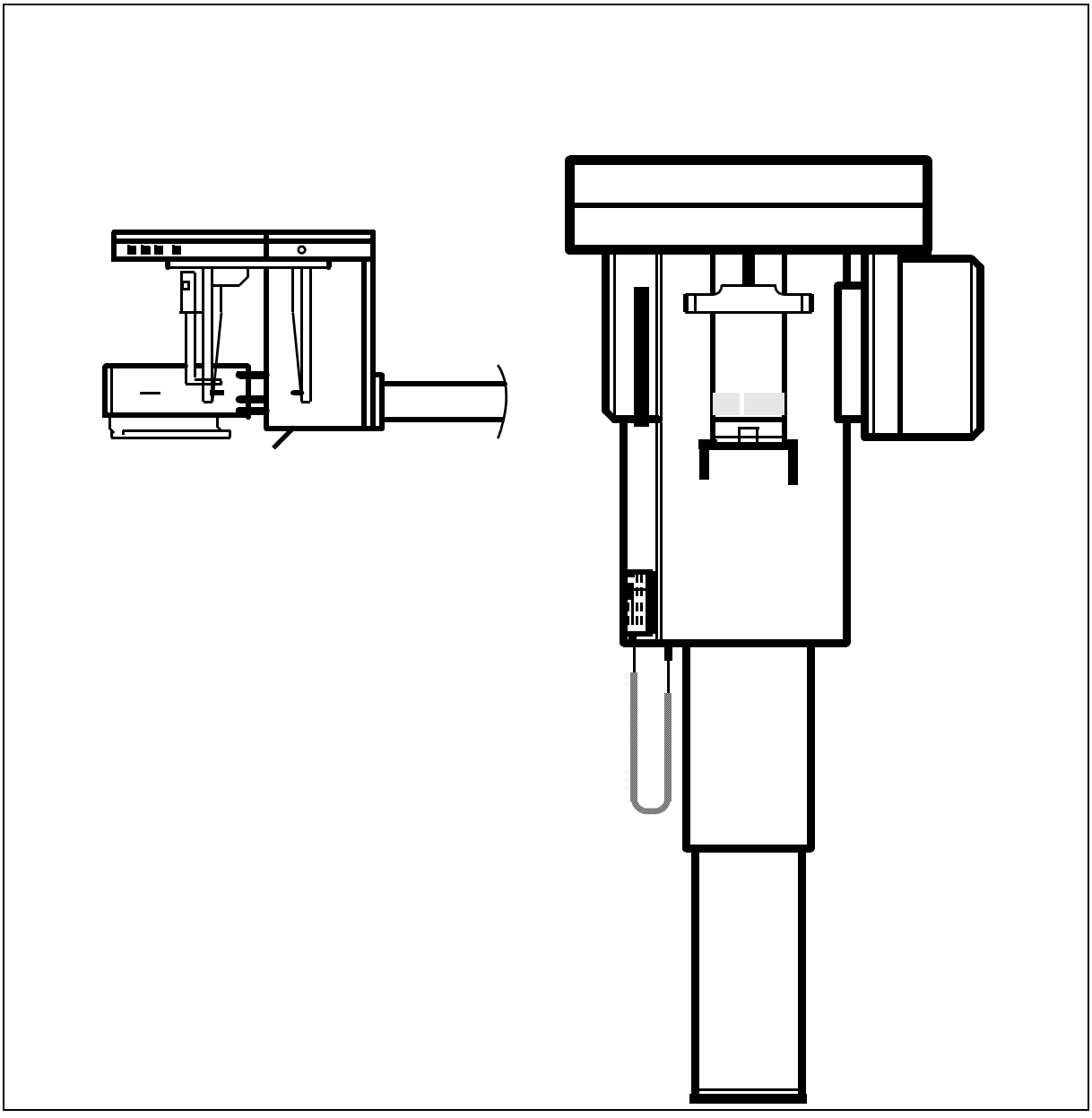ATTENTION

Interference with electromedical devices by radio telephones:
To guarantee the operational safety of electromedical devices, the operation of mobile radio telephones in the medical practice or hospital area is prohibited.

ATTENTION

When opening the equipment:
Please observe the safety measures for handling PC boards.
Touch a ground point to remove any personal electrostatic charge before touching the components.

New since: 11.2004

Modification compared with last edition: 10.1999

Chapter or paragraph, page
5 “Concluding work” on page 9

Contents

1 Remove covers 3


2 Replace software 4

3 Replace the primary diaphragm 5

4 Select service routine S.21 6
4.1 Prepare service routine 6
4.2 Programming the diaphragm number in service routine 21 7

5 Concluding work 9

1 Remove covers

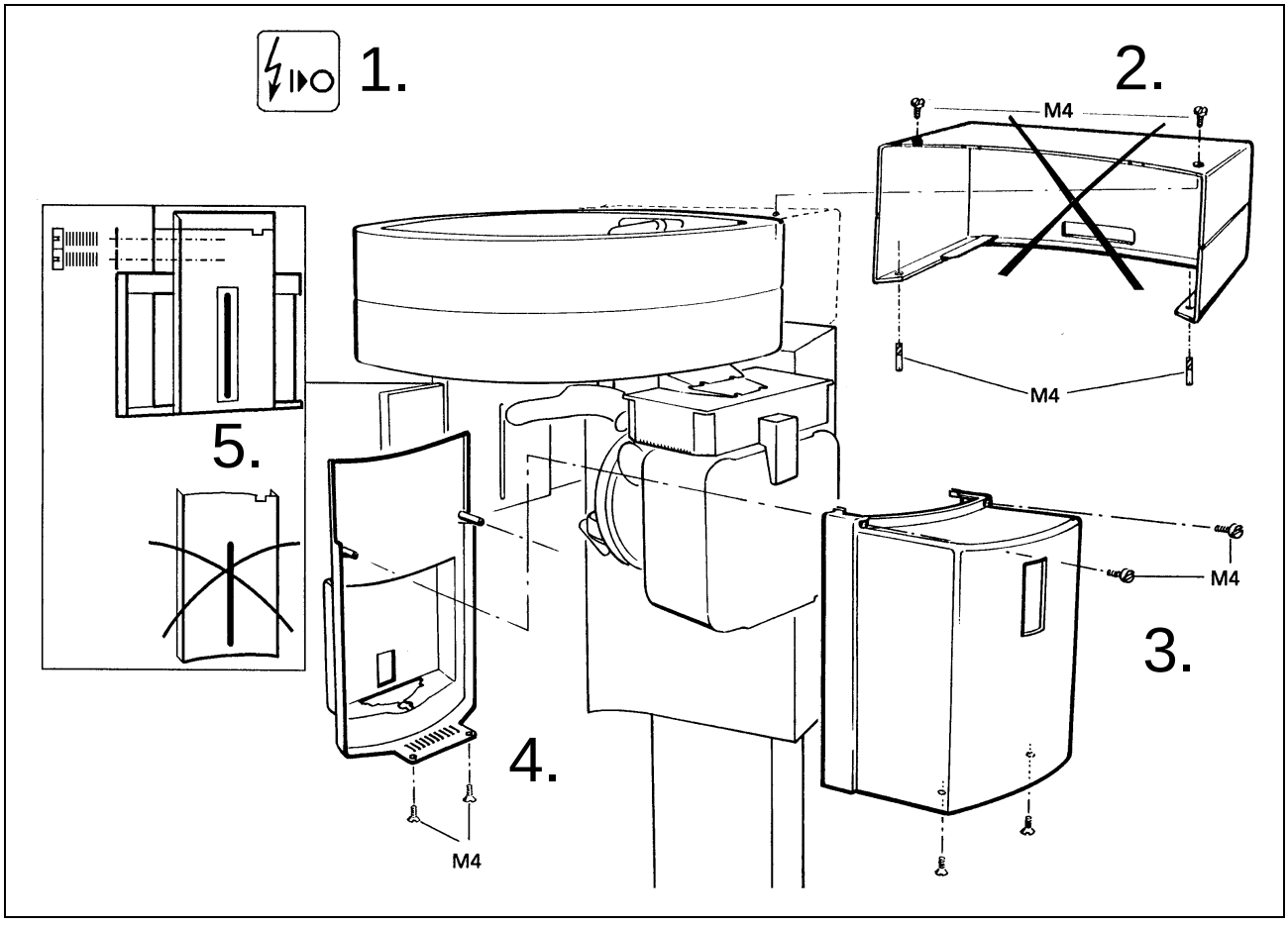
NOTE

The* Service Manual *is an essential requirement for the installation.

1. Disconnect the unit from the power supply.
2. Unscrew two **M4 setscrews** from below.
 Unscrew two **M4 cheese head screws**.
 Remove the upper cover (no longer required).
3. Unscrew two **M4 cheese head screws M4**.
4. Unscrew two **M4 countersunk head screws** and remove the front cover.
5. Unscrew the scattered radiation protection (no longer required).
 Fasten the supplied scattered radiation protection.

English

2 Replace software

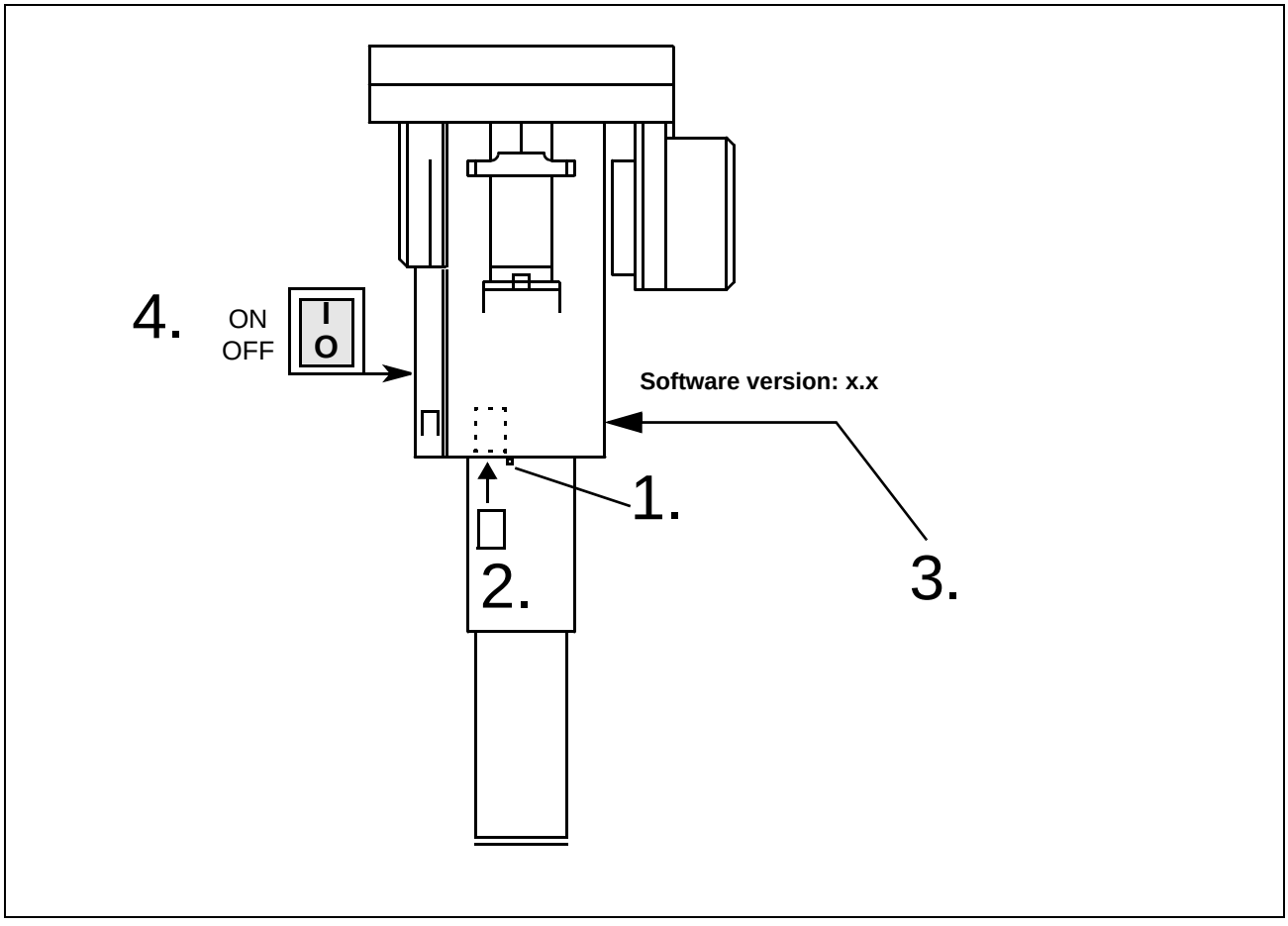
NOTE

The unit must be switched off.

1. Push the ejection button of the PC card and remove the card.
2. Insert the supplied PC card with new software (unlabelled side must point to the unit).
3. Stick new label “Software version x.x” over the old label.
4. Switch on the unit.

- The display shows **E03.52**.
- Verify the unit identification with Service Routine **S.17** (see Service Manual):
 - Press the Memory key
 - Press the Return key **R**
 - Briefly press the Service key

3 Replace the primary diaphragm

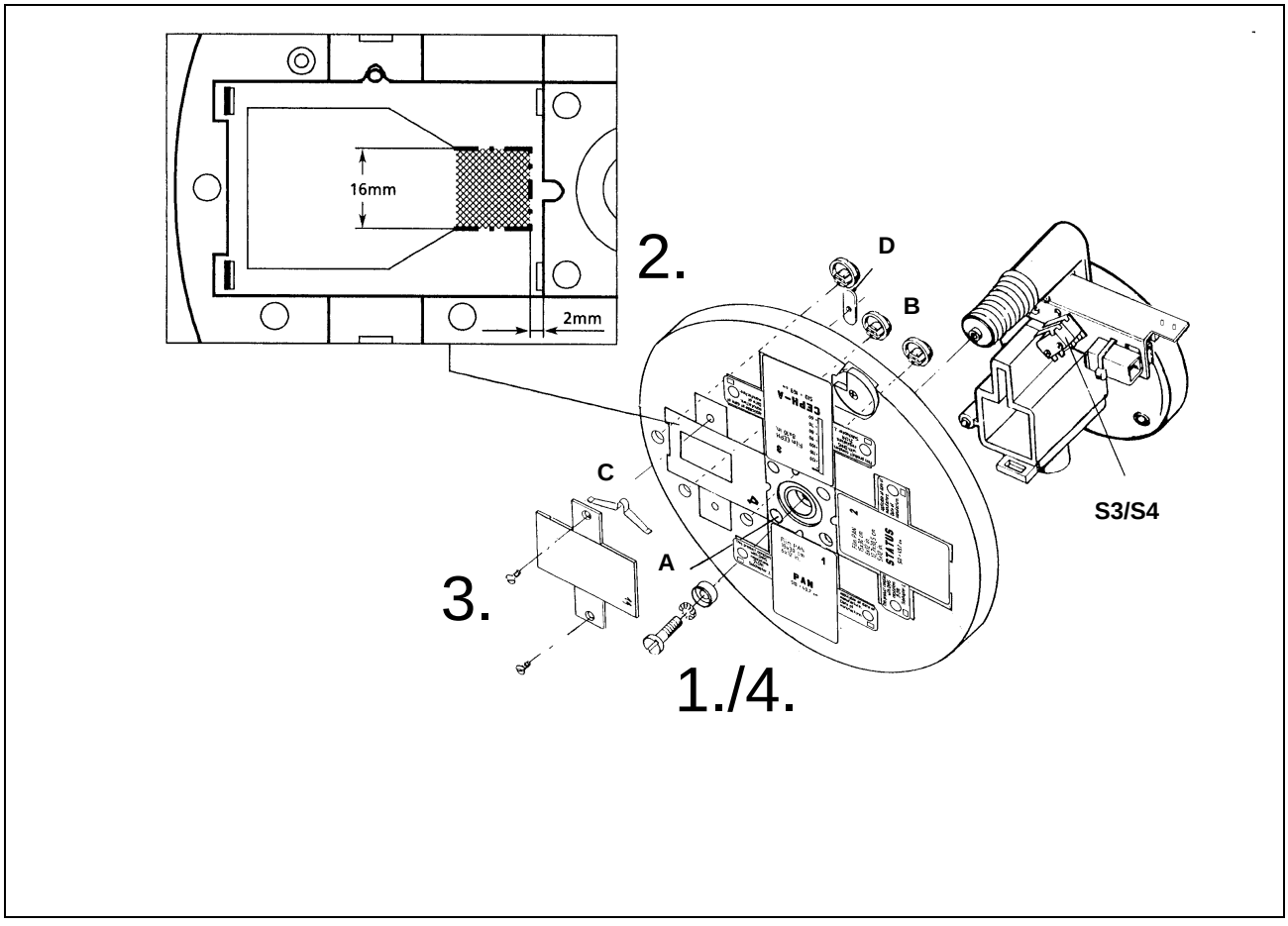1. Unscrew the diaphragm wheel (countersunk head screw, serrated washer and bush).

- If all diaphragm locations are occupied, unscrew the diaphragm 2, 3 or 4 according to the customer's wishes.

2. File down the hatched surface on the diaphragm wheel.

3. Insert three eccentrics **B** (if not present).
Insert spring **C** (if not present).
Fasten supplied primary diaphragm **11** with two threaded plates **D** and countersunk head screws.


4. Reattach diaphragm wheel with countersunk head screw, serrated washer and bush.

ATTENTION
When fitting the diaphragm wheel take care that the switch S3/S4 is not damaged by the coding pins.

- The supplied unit passports and envelopes are required only in the USA / CANADA.

English

4 Select service routine S.21

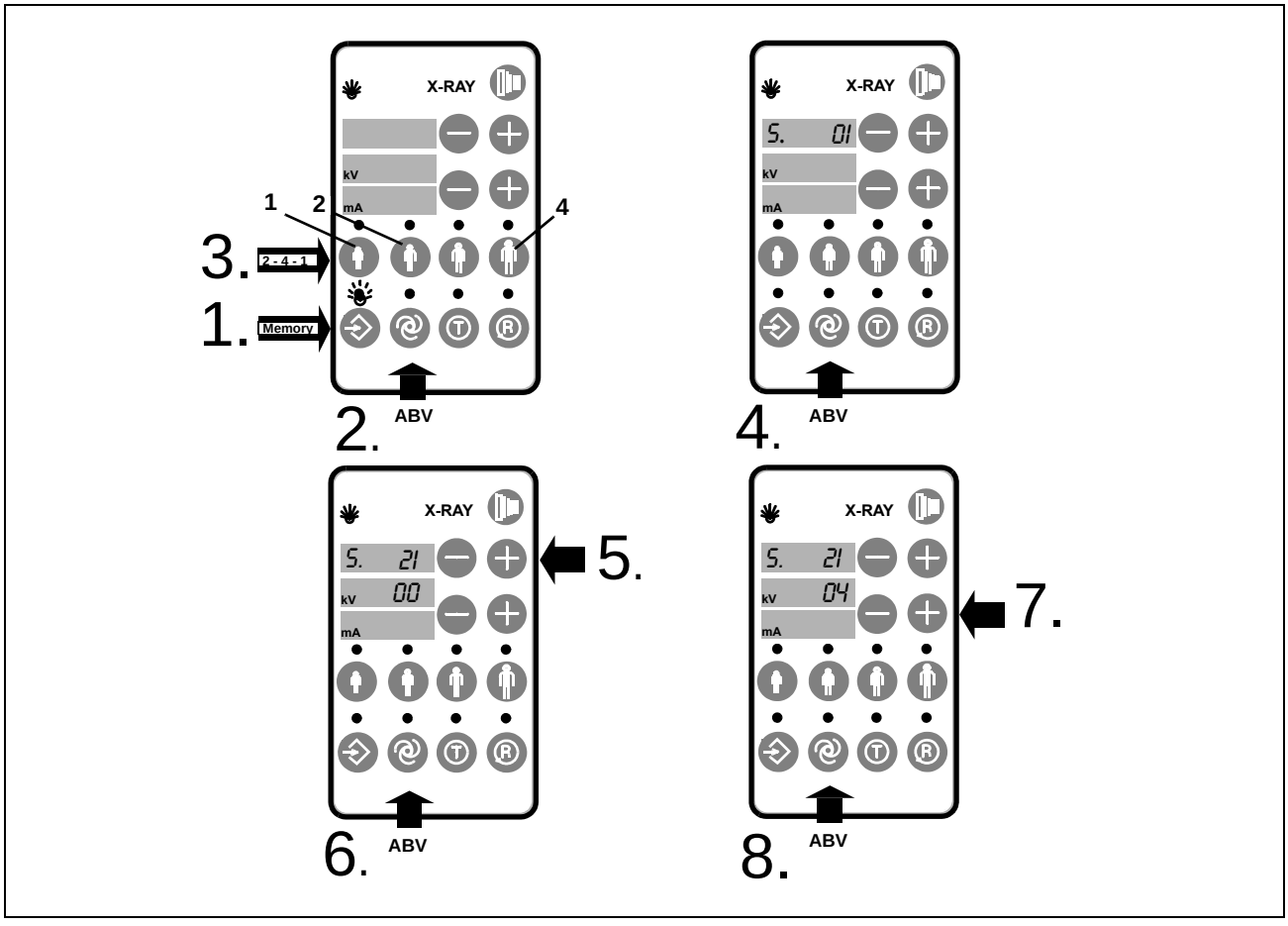4.1 Prepare service routine

- Perform the service routine **S.21** after installing the diaphragm.

1. Press the **Memory** button (LED above the button flashes).

2. Then press the **ABV button** (approx. 4s) until the digital displays disappear.

3. **Within 3s**:
Press the **patient symbol buttons** in the sequence **2 – 4 – 1**. The service mode is now selected.

- If you do not adhere to the sequence when selecting the service routine, press a wrong button or exceed the time (see 3.). The system then switches automatically to normal mode.


4. Press the **ABV button**.
Service routine **S.0 1** appears on the upper digital display (program display).

5. Select the service routine **S.21** with the **“+” (program selection)** button.

6. Press the **ABV button**.
00 appears on the kV display

7. Press the **“+” (kV/mA)** button several times until the service code **21** appears on the kV display.
This step must be performed for safety reasons in all service routines which can give rise to a hazard.

8. Press the **ABV button once again**. The digital displays flash.
The diaphragm location is displayed on the kV display (01 to 04). If the diaphragm wheel has not engaged, then 00 appears on the kV/mA display.

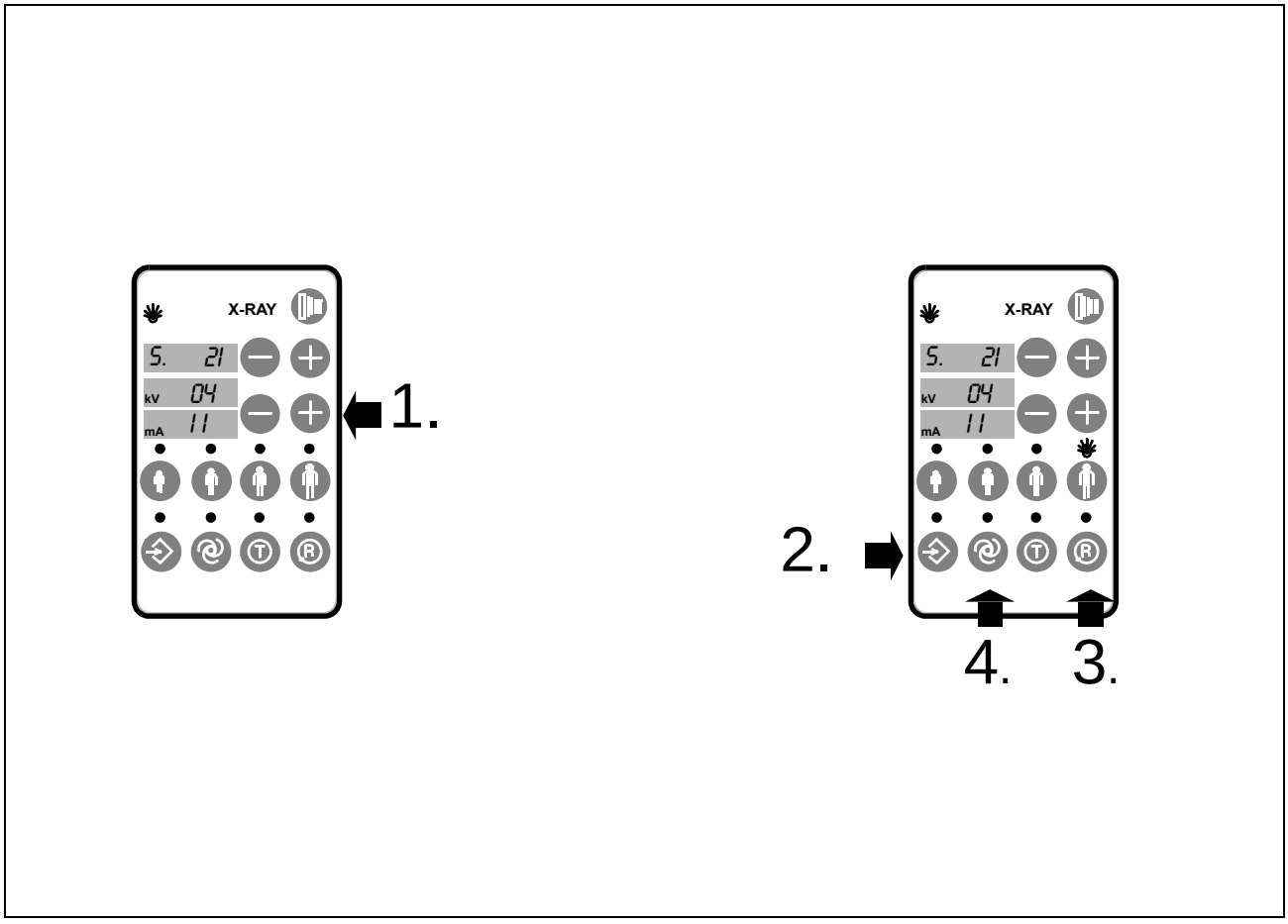English

4.2 Programming the diaphragm number in service routine 21

- The transverse diaphragm must be in the radiation position, since otherwise no change is possible.


NOTE

In the following collimator position 4 is described by way of example. According to the position of the diaphragm to be replaced, other positions can also be referred to.

Case 1: Previously there was no diaphragm in the location

- Display for diaphragm location 4 (example):
 kV 04
 mA 00

1. Enter the value **11** in the mA display with th **"+" (kV/mA)** button.
2. Press the Memory button.
3. Press the Return button **R**.
- The patient symbol LED of the corresponding collimator position lights up.
4. Press the **ABV button** briefly. The service routine S.21 can be exited.
- **To exit the service routine switch the unit OFF and back ON.**

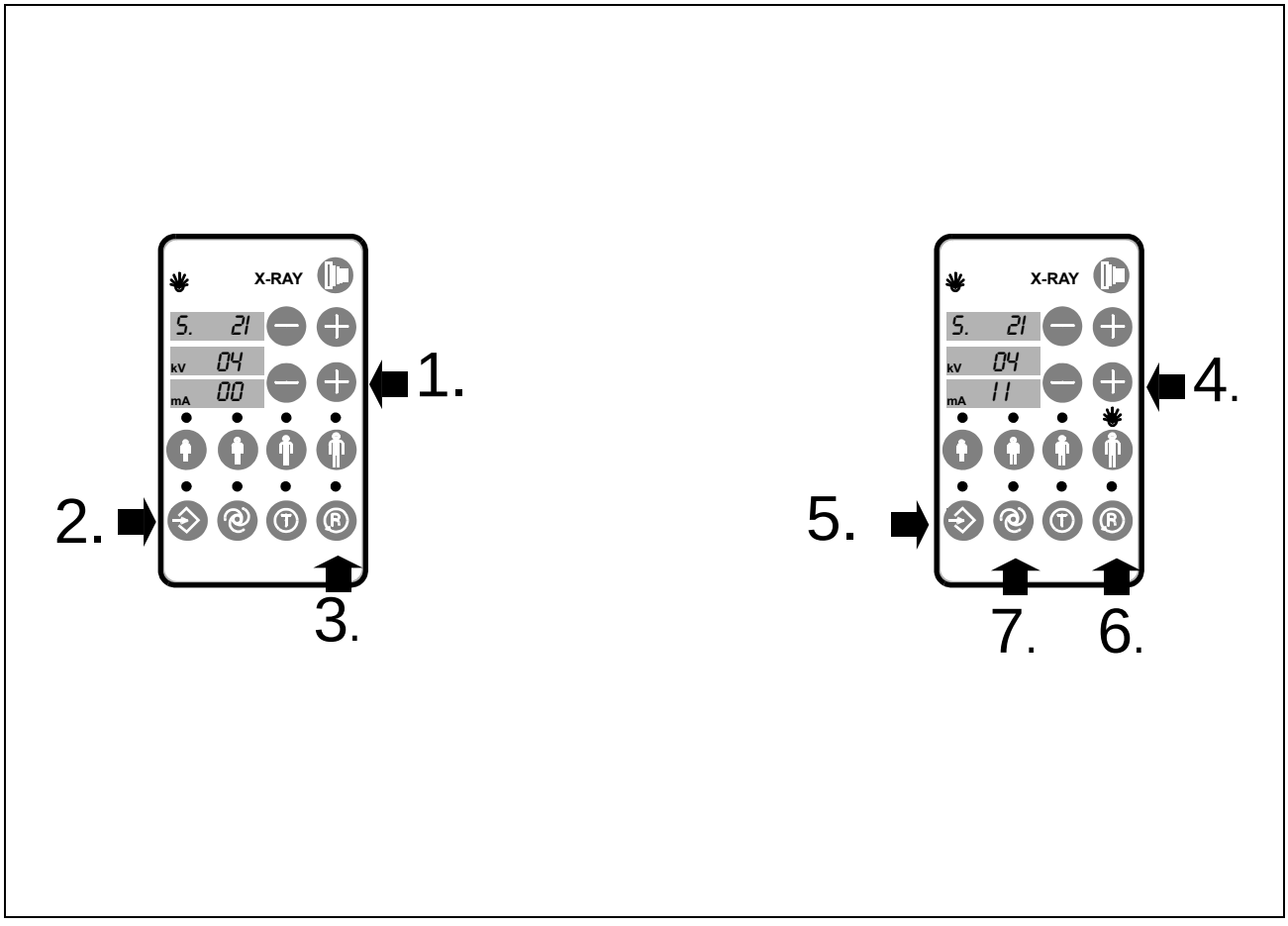Case 2: A diaphragm was previously already present in the location

- Display for diaphragm location 4:
 kV 04
 mA xx

1. Enter the value **00** in the mA display with th **“+” (kV/mA)** button.

2. Press the Memory button.

3. Press the Return button **R**.

4. Enter the value **11** in the mA display with the **“+” (kV/mA)** button.

5. Press the Memory button.

6. Press the Return button **R**.

- The patient symbol LED of the corresponding collimator position lights up.

7. Press the **ABV button** briefly. The service routine S.21 can be exited.

- **To exit the service routine switch the unit OFF and back ON.**

5 Concluding work

- Attach the sheet metal cover and paneling **A**, referring to the **Service Manual**.
- Check and adjust the diaphragm **11** as described for diaphragm 1 in the **Service Manual**, ”**Check and adjust the X-ray beam for panoramic exposures**”.
- Attach the covers and use the supplied upper cover.
- Insert the diaphragm insert in the film holder and swing diaphragm 1 into place.
- Enter the new diaphragm in the unit passport.
- Perform test runs with the TSA programs. Make sure that the procedures run noise-free and collision-free (actuator sleeve and cover).
- Hand the supplied accessories and the operating instructions over to the customer.
- **The unit is now ready.**

We reserve the right to make any alterations which may be required due to technical improvements.

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